


HOT DOG AIR-FREE PATIENT WARMING



**MAINTAIN NORMOTHERMIA
WITH THE TOUCH OF A BUTTON**

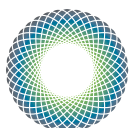

**LATERAL
MEDICAL**



WC71 - Single Port Temperature Management Controller

- 21cm W x 29cm H x 18cm D
- 2.9kg
- 17.5cm Touchscreen with Intuitive User Interface
- Controls **one** Warming Blanket
- Temperature Monitoring*
- Integrated Slide Shows
- Positioning Guides
- Includes IV Pole Clamp
- VESA Mount Compatible (75mm)
- Integrated Self Testing



WC77 - Multiport Temperature Management Controller

- 21cm W x 29cm H x 18cm D
- 3.2kg
- 17.5cm Touchscreen with Intuitive User Interface
- Controls up to **four** Warming Blankets and **one** Warming Mattress
- Temperature Monitoring*
- Integrated Slide Shows
- Positioning Guides
- Includes IV Pole Clamp
- VESA Mount Compatible
- Integrated Self Testing

The WC77 Multiport Controller provides maximum versatility and efficiency for patient temperature management.


WC77 Back



Controller Features


Integrated Training Slideshows

Controllers include educational slideshows to help you get started and learn the BEST Practices.

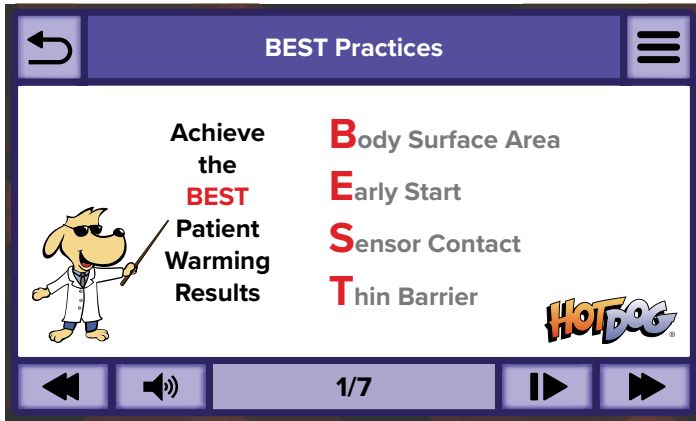


Temperature Monitoring & Auto Mode*

Temperature monitoring uses patient temperature to automatically control the warming devices.

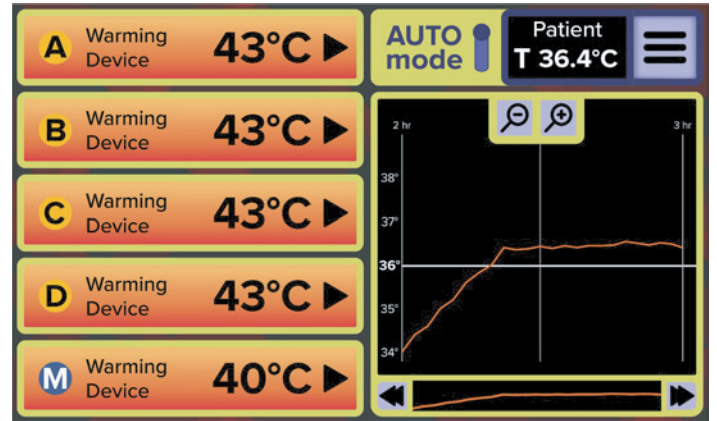

Advanced Temperature Control

Control all warming devices individually in 1°C increments or together with the touchscreen.

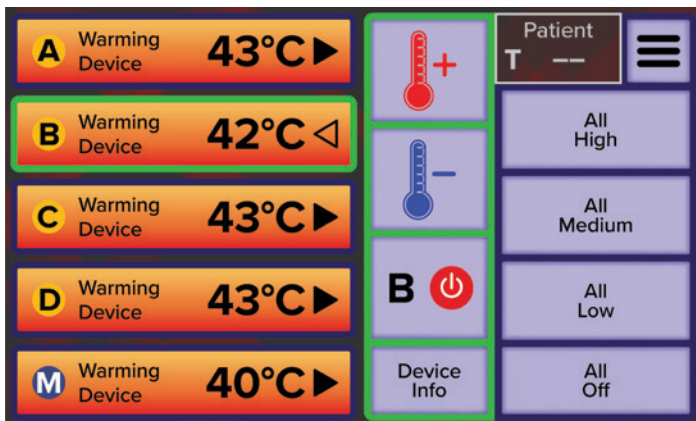

Positioning Guides

Includes visual blanket guides for various surgical positions to make set-up easy for new users.



The "smart" controllers deliver low-voltage electricity d by an advanced sensor system. There is a 3-year warranty for controllers.

*Temperature Monitoring feature is CE-marked and pending FDA clearance.

Warming Blankets


B500 Universal Warming Blanket

43.2cm W x 76.2cm L
0.8kg



HotDog warming blankets are lightweight, flexible, and designed to meet all of your warming needs. The outer shell is a soft, non-porous urethane embedded with an anti-microbial agent that prevents the growth of micro-organisms — designed for easy cleaning. The semi-conductive fabric produces uniform heat output over the entire surface of the blanket. Temperature is controlled from 37-43°C. The blankets have a 2-year warranty and a 30-month expiry.

HotDog can help maintain normothermia for patients of all ages. See page 6 for our paediatric product offering.

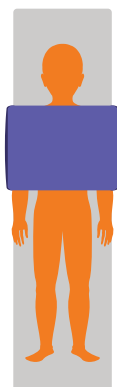

1 B500




1 B500



2 B500s



2 B500s
(1 Folded)

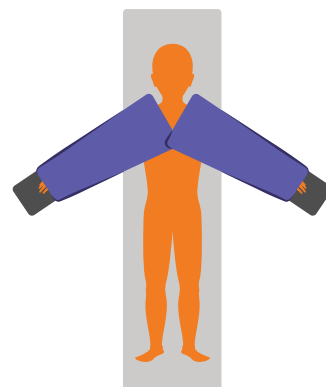

2 B500s

The B500 is the most advanced and versatile warming blanket. It has rounded edges for a softer touch and formable edges to maintain shape. It also includes a more flexible heater fabric for better drapability and effectiveness.

Warming Blankets

B103 Lower Body Warming Blanket

74.9cm W x 98.4cm L
1.36kg




B104 Full Body Warming Blanket

74.3cm W x 128.6cm L
1.81kg




B105 Multi-Position Warming Blanket (Two tethered panels)

50.8cm W x 90.63cm L (x2 Panels)
1.45kg

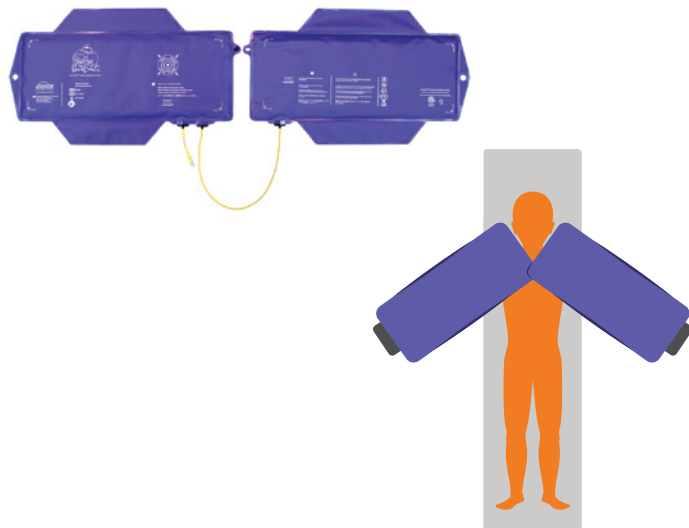

B107 Head Warming Blanket

19.7cm W x 70.5cm L
0.54kg

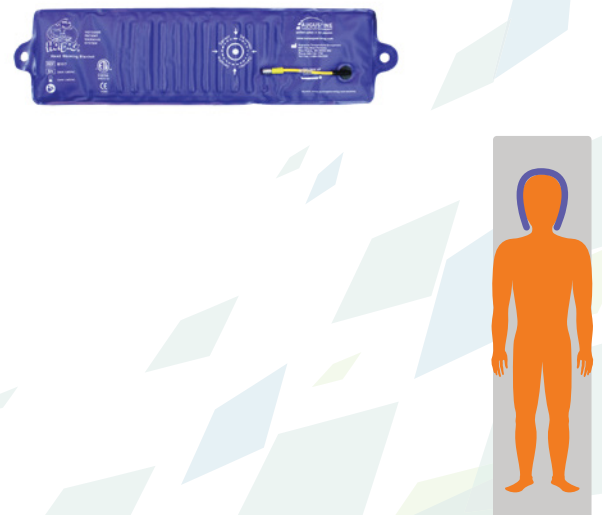

Warming Mattresses

U101 Underbody Warming Mattress

82cm

49.53cm W x 81.28cm L

0.9kg



U102 Underbody Warming Mattress

127cm

49.53cm W x 127cm L

1.5kg



U300 Trendelenburg Position Warming Mattress

58.42cm W x 99.7cm L

1.04kg



HotDog warming mattress overlays are designed to provide safe warming under the patient (36-40°C settings). The fabric heater is layered over pressure-relief foam, providing a uniformly heated, pressure-neutral pad between the patient and the table surface. The outer shell is a non-porous urethane embedded with an anti-microbial agent that prevents the growth of micro-organisms — designed for easy cleaning. The mattresses have a 2-year warranty and a 30-month expiry.

Mattresses are mostly radio-lucent. Only the bus bars along the long edge, the sensor, and connecting wires are radio-opaque.


Paediatric Warming

B270 Small Paediatric Head Warming Wrap

For patients <8kg

11.43cm W x 33.02cm L

0.09kg

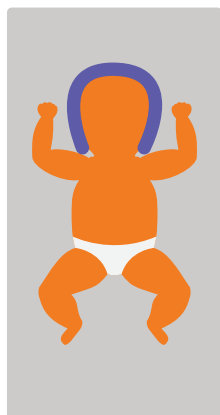


B271 Large Paediatric Head Warming Wrap

For patients >8kg

13.34cm W x 38.74cm L


0.18kg

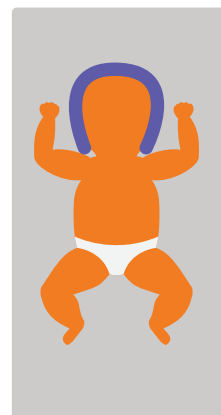

B203 Paediatric Lower Body Warming Blanket

45.72cm W x 55.88cm L

0.45kg

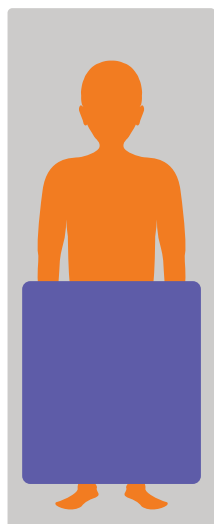

U220 Paediatric Underbody Warming Mattress

For patients 20kg

44.45cm W x 73.66cm L


1.0kg



WaffleGrip® System

WaffleGrip Trendelenburg Positioning System

- Technology integration that optimises surgical efficiency, safety, & outcomes
- Increase warming surface area by 3.5x for improved patient outcomes
- Secure the patient with WaffleGrip's unique anti-slip technology
- Incorporated slide sheet allows effortless positioning to improve staff safety
- Perineal drape helps manage fluids

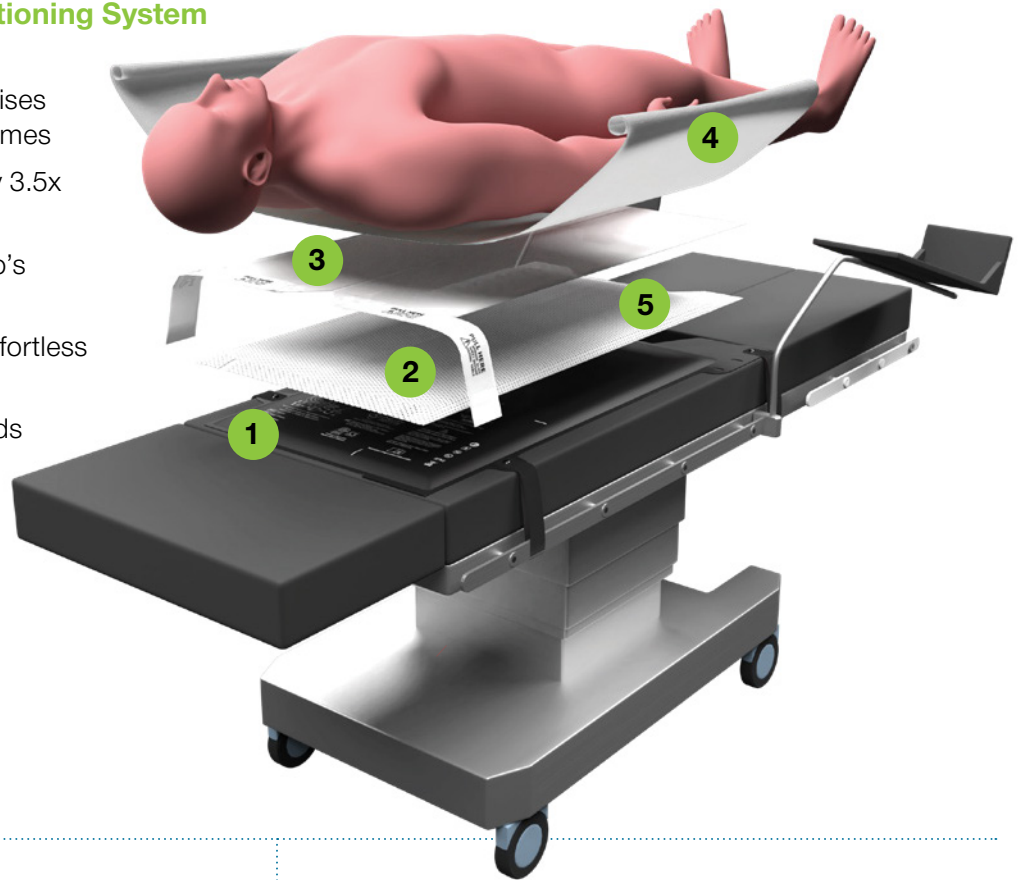

U300 Trendelenburg Position Warming Mattress

58.42cm W x 99.7cm L
1.04kg



WaffleGrip must be used with the U300.

A301 WaffleGrip Trendelenburg Positioning System + BackSaver Slide Sheet, 10 Pack

Each A301 Kit Includes:

- 1 WaffleGrip Pad
- 2 BackSaver Slide Sheet
- 3 Drawsheet with WaffleGrip
- 4 Perineal Drape (Hidden)
- 5 Disposable Strap (not pictured)


WaffleGrip
Warm & Secure

Accessories

A101
Blanket Cable
4.0m



A112
Mattress Cable
5.0m



A116
Roll Stand
(Includes 2 wire storage racks)




A108
Wire Storage Rack
IV Pole Mount



A109
Wire Storage Rack
Wall Mount



A110
Wire Storage Rack
WC7X Clamp Mount
(Only for use with WC7X)

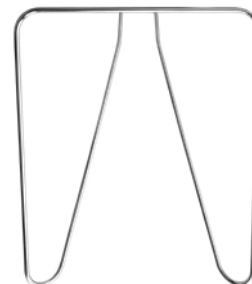

CV10
Clinician Warming
Vest Heater



C21
Clinician Vest
(XS-S-M)

C22
Clinician Vest
(L-XL-XXL)



C301
Clinician Vest Kit
(XS-S-M)
Includes C21 Vest, CV10 Heater
C3021
Clinician Vest Kit
(L-XL-XXL)
Includes C22 Vest, CV10 Heater



HotDog patient warming is designed to prevent and treat hypothermia in surgical patients. It can be used before, during, and after surgery. The complications associated with unintended hypothermia include increased rates of: bleeding, surgical site infections, adverse cardiac events, and morbidity. Rx Only.

Regulatory & Quality Standards Met

FDA 510K clearance for adult and pediatric use
(K052392, K092807, K094056, K112488, K201779)
IEC 80601-2-35 Medical Use Heating Pads
IEC 60601-1 Class II Type BF
IEC 60601-1-2 Emissions and Radiation CE Mark
ISO 13485 registered
ETL listed (US & Canada)

Patents:

Augustine Surgical Inc. (ASI) and related parties have been issued sixteen patents by the United States Patent Office containing approximately 150 claims. Twenty-two other patent applications are pending. One PCT patent has been issued. Five additional patents have been issued for securement solutions. The patents and applications, ASI believes, cover every safe and effective way to utilise flexible heating fabric in a surgical setting. For more information on patents visit hotdogwarming.com/patents.

A - Air Free

B - Better Warming

C - Cost Effective



LATERAL
MEDICAL



Augustine
Surgical Inc.
Making Surgery Safer

Contact: 1300 144 734 or
sales@lateralmedical.com